

## How to use Colgate Water Flosser



Plug in the nozzle tightly until you hear a click



Open the lid, pull the water tank down until you hear a click & fill in with water.



Extend nozzle & glide along teeth

Do not point on gums



Press on/off button\* and select your preferred mode

Default mode is gentle mode

\*Do not aim and spray water directly under the tongue, into eyes, ears, nose, or other parts of the body. Pressure from the water stream may cause damage.

## Nozzle assembling and disassembling



Assemble the nozzle into the water flosser nozzle jack. Insert firmly until "click" sound is heard.



Press and hold the nozzle button for disassembling. Gently pull out the nozzle.



To store, place the nozzle in the storage compartment.



**REMINDER:**  
Replace your water flosser nozzle every 6 months



Retain your receipt for warranty. Contact Consumer Affairs for any product queries & assistance. Toll-free number: 1-800-883-918

Colgate®

NEW

Ask your dentist TODAY!

# WATER FLOSSER

Floss your way to a better smile!



# Floss your way to a better SMILE!



**Blast Away Plaque :**  
Between teeth  
and along the gums



**Up to 60 seconds of  
Water Capacity :**  
Easy-to-fill  
water tank




**3 Pressure Modes :**  
Personalize your setting  
Gentle, Regular, Pulse



**Rechargeable Battery :**  
More than 4 weeks' battery  
life with daily use



**Mind The Gaps :**  
Easy Clean  
around braces

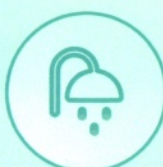


**Waterproof Design**  
Suitable for usage  
in the shower

## Ultra fine stream of water to clean between teeth



**Personalize your  
setting with 3 Modes**

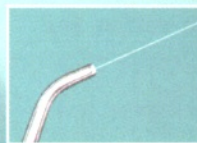

**Gentle mode**  
Light and soft water pressure  
Friendly to new users



**Regular mode**  
Thorough wash  
Effective cleaning



**Pulse mode**  
Targeted stream of water  
released at one time

