


Prontosan[®] Wound Irrigation Solution

The solution for acute, chronic and burn wounds

B. Braun Wound Care. Right. From the start.

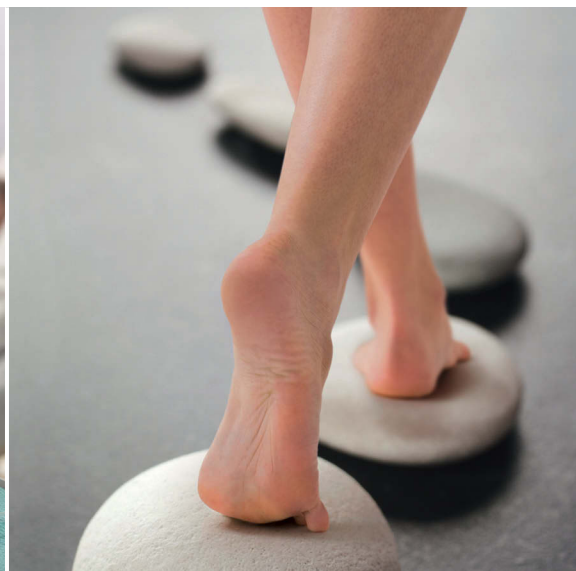

Wound Management

B | BRAUN
SHARING EXPERTISE

Prontosan® Wound Irrigation Solution

The solution for acute, chronic and burn wounds



Removes & prevents biofilm

Prevents infections

Reduces healing time

Painless dressing changes

Well-known substances with low allergenic potential

Compatible with commonly used wound dressings

Can be used up to 8 weeks after first opening

Prontosan® contains betaine and polihexanide.

Betaine reduces the environmental surface tension of the biofilm which results in a maceration of the biofilm and finally allows a better penetration into it. The mechanism of action can basically be compared with a dishwashing detergent.

Polihexanide is an antimicrobial agent with broad spectrum biocidal activity. The mode of action can be described as a non-specific electrostatic interaction with the bacterial cell wall. The attachment of polihexanide to the bacterial cell wall results in a disorganisation of the biological structure of the bacteria.




USE OF PRONTOSAN® IN NPWT

The negative pressure wound therapy (NPWT) irrigating systems provide an even more effective means for delivering solutions such as Prontosan® Wound Irrigation Solution to complicated wounds with the opportunity for:

- › Preventing infection
- › Decreasing healing time
- › Reducing the risk of amputation

CAUTION: Use on hyaline cartilage and in aseptic joint operations is contraindicated.

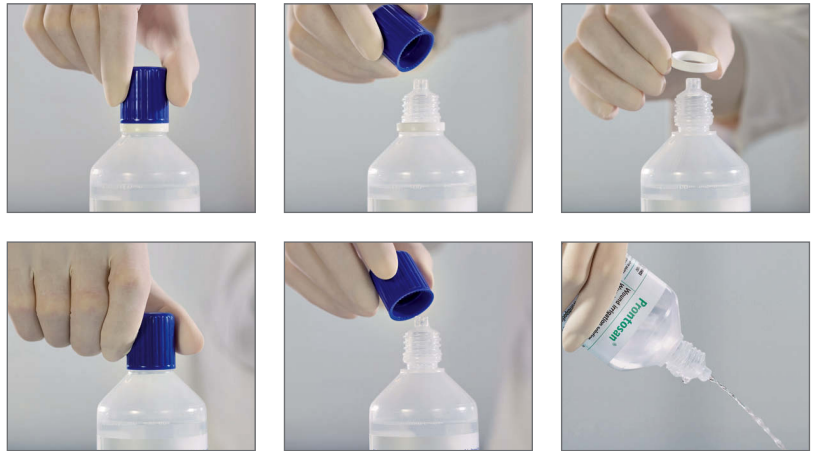

Indications

Prontosan® Wound Irrigation Solution can be used for the prevention of an infection in acute and chronic wounds including

- › Traumatic wounds
- › Postoperative wounds
- › Chronic skin ulcers (e.g. venous, diabetic or pressure ulcers)
- › Burns 1st and 2nd degree

Before Application

- › To open the bottle remove distance ring and screw the cap clockwise

Advantages

- › Ready to use
- › Easy and simple to use
- › Direct application from the practical squeeze bottle
- › Well-known substances with low allergenic potential
- › Can be warmed up to body temperature
- › Can be used on sensitive or irritated skin
- › No risk of further contamination of wounds
- › Can be used in combination with standard and advanced wound dressings



Prontosan®	Size	Pcs/Pack	Article number
Pod	40 ml	24	400484
Bottle	350 ml	1	400403
Bottle	1,000 ml	1	400446

Composition: 0.1% Polyaminopropyl Biguanide (Polihexanide)
Betaine Surfactant

