IV Rail TOSO

New Delac Parts for IV

Product Name		Product Code
New Delac IV Roller Allowable load: 10 kg (22 lb)		00210052
New Delac End Stop B	White	00210332
	Natural	00210135
New Delac Splice (for Alumi Rail)		00210128

Supertrac Parts for IV

Product Name	Product Code	
Supertrac IV Roller For a curve rail, use it too ther with a rail of 300 mm (12"R) or longer.	00240001	
Supertrac End Cap	00240073	
	Natural	00240019

Product Name	Product Code	
Supertrac End Stop	Natural	00240020
Supertrac Splice (for Alumi Rail)		00240004

IV Hooks

Product Name		Product Code
IV Bottles Hook Allowable load: 10 kg (22 lb)	F	01080010
IV Bottles Hook 30GN Allowable load: 30 kg (66 lb)		01080012
IV Wall Hook Allowable load: 10 kg (22 lb)		01080013

Product Name	Product Name	
IV Bottles Hook	A600-P1	0045000
	A800-P1	0045000
IV Bottles Hook	A600-P3	0045001
	A800-P3	0045001
IV Bottles Hook	B1000-P5	0045001
	B1200-P5	0045002

201001E




IV Rail




- •IV rail designed for easy use and safety.
- Ceiling-suspended for easy attachment and removal of IV bottle hangers.
- •Various combinations possible for a wide range of applications.

Conventional IV Rail

(79")

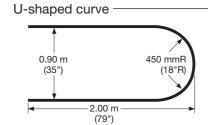

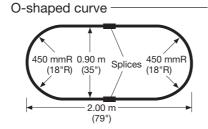

New Delac (for hospital ward applications)



IV Roller Dimensions




Materials and Finishes

Parts		Material	Finish
Rail Aluminum alloy extrusion material		,	Alumite
IV Roller Wheel		Acetal resin	
iv Holler	Hook	SUS 304 5 Ø	

New Delac Rails

22.5 mm (.89") 1 20 3 1 15.8 mm	22.5 n
106 mm (4.2") 90 mm (3.5")	
21 777	——→ 7.0 m (.28'
31 mm (1.2")	

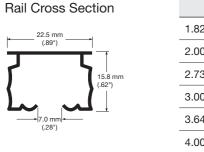


	Size	Product Code	
	SiZe		Natural
1.82 m (72")		00210346	00210120
2.00 m (79")		00210347	00210121
2.73 m (107")		00210348	00210122
3.00 m (118")		00210349	00210123
3.64 m (143")		00210350	00210124
4.00 m (157")		00210351	00210125
Curved Rail (180 mmR/ 7.1"R)	0.40 m x 0.40 m (16" x 16")		00210126
	1.00 m x 1.00 m (39" x 39")		00210035


Supertrac (for ICU and CCU and operating room applications)





IV Roller Dimensions Rail Cross Section

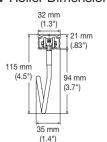

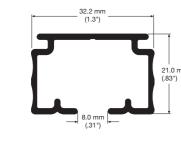

Materials and Finishes


Parts		ts	Material	Finish	
	Rail Aluminum alloy extrusion material		Alumite		
	IV Roller	Wheel	Special resin		
	IV Roller	Hook	SUS 304 6ø		

Supertrac Rails

Size	Product Code	
Size	White	Natural
3.00 m (118")	00240087	00240023
4.00 m (157")	00240088	00240026
Curved Rail (500 mmR/20"R) 1.00 m x 1.00 m (39" x 39") (Ordered separately) *1		00240058

^{*1} Curved 1.00 m x 1.00 m (39" x 39") rails (ordered separately) are 500 mmR (20"R). Note: Possible radius sizes····· Natural: 450 mmR (18"R), 800 mmR (31") or more

Common Parts

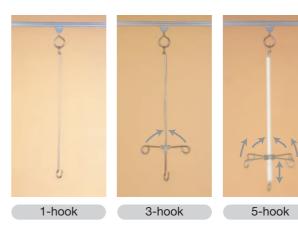

IV Bottles Hooks

1-hook, 3-hook and 5-hook hangers are available, each in two lengths, for a total of six types. The hanger length of the 5-hook type can be adjusted. And the hooks on the 3- and 5-hook types are foldable for convenience when carrying and storing.

Specifications

Products	Hooks	Extension Range	Length	Allowable Load	Materials and Finishes
B1000-P5	5-hook	638 mm (25") \iff 1000 mm (39")	1000 mm (39")	Lower hook: 5 kg (11 lb)	Aluminum, Acetal,
B1200-P5	3-1100K	738 mm (29") 😂 1200 mm (47")	1200 mm (47")	Intermediate hooks: 1kg (2.2 lb)	Zinc die casting
A600-P3	3-hook		600 mm (24")	Lower hook: 5 kg (11 lb)	SUS 304 6ø
A800-P3	3-1100K		800 mm (31")	Intermediate hooks: 1kg (2.2 lb)	
A600-P1	1-hook		600 mm (24")	Lower hook: 5 kg (11 lb)	505 304 60
A800-P1	1-1100K		800 mm (31")		


IV Wall Hooks

Slow Return function slowly retracts the hook.

For easier removal of intravenous bottles or bags, the hook retracts slowly so that it does not impede the hand.





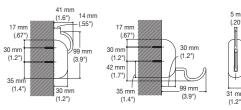

Precautions when mounting IV Hooks

*Securely attach each hook with two of the supplied 3.5 mm (.14") x 25 mm (.98") pan-head screws.

*To assure hooks remain firmly attached, mount on plywood at least 12 mm (.47") thick or on another substrate of the same

Precautions Concerning IV Hooks

Allowable load: 10 kg (22 lb)

*Securely attach each hook with three of the supplied 4.0 mm (.16") x 40 mm (16") pan-head screws.

The weight of an IV bottle causes the runner to lock, but it unlocks if the IV hanger is gently lifted, allowing the runner to move smoothly. Special plastic used for the runner wheels assures the runners are extremely quiet.






Locked

IV Hook Dimensions

<Fixed type> <30GN fixed type>

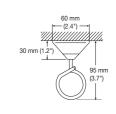



30GN fixed type features

- *As a safety feature, the hook comes apart for when loads of 30kg (66lb) or more are applied.
- *Be aware that even with loads lighter than 30kg, dynamic load may cause a hook to become detached.

 $\mathbf{3}$